


LAF Software Update and new Configuration

Step1

- You need a USB stick.
- **Very important** the USB stick must be called **UPDATESTICK** written together, and all letters capitalized.
- they do their new software or configuration.
on the USB stick make sure the file is not in a folder

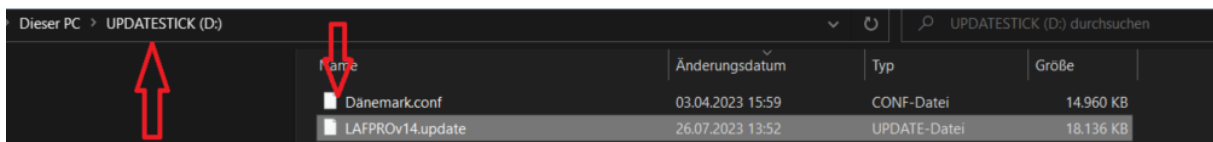


Step2

- insert the USB stick into the port.
- on the left is the bar and on the right is the Classic displayed




Step3

- Go to the technical settings.
- Press Software Update

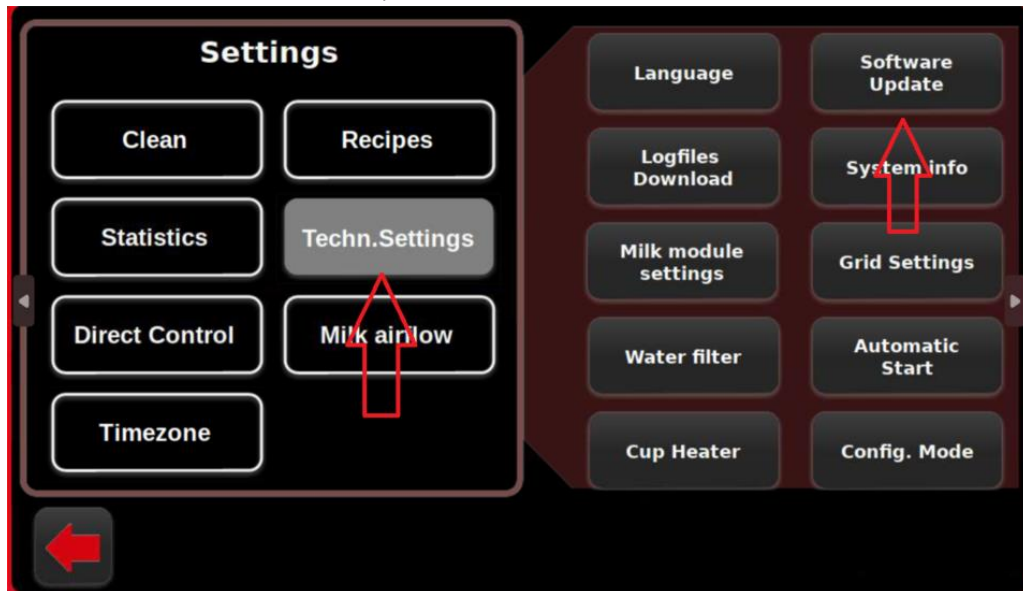

Step4

- first insert the USB stick and then press Update.

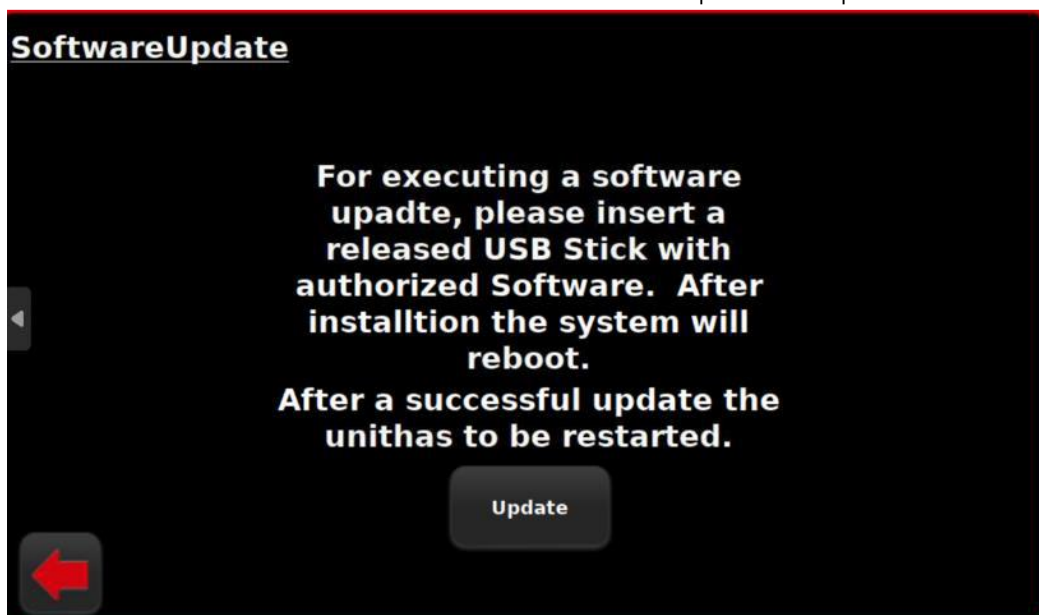

Step5

- now a white screen should appear with the respective files that are on the USB.
- select the respective file and then press Install.
- as soon as the file has been successfully uploaded a green check mark appears.
- now you can turn off the machine and then restart it the new update or configuration should be installed on the machine now.

